Body TRAINER

- Elettrostimolatore cutaneo • Cutaneous electrostimulator
• Électrostimulateur cutané • Elektrischer Hautstimulato
• Electroestimulador cutáneo

MADE IN ITALY

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Cutaneous electrostimulator

Cutaneous electrostimulator

- 1 180° application head**
- 2 Ergonomic handle for all hands**
- 3 Intensity control wheel on the body**
- 4 Indicator window showing the intensity level (1-10)**
- 5 Led light indicating the device functioning**
- 6 ABS resistant impact body**
- 7 Battery compartment cover (1 x AA 1,5 V)**
- 8 Protective cloth case**

Please carefully read all the instructions before use.

Cutaneous electrostimulator

We all want to be always in a good shape, all days; the light malaises of everyday most of the time depend on our muscle tension, on minor soreness and on everyday stress.

Body Trainer is a portable non-invasive electronic device, to relieve the pain naturally without using drugs. It is effective against muscle contractions, to relax muscles and against the everyday stress.

Its efficacy is based on the medical technique called "Transcutaneous Electrical Nerve Stimulation" (TENS) that through an electronic circuit transmits low potential electrical impulses that stimulate some particular skin corpuscles (skin receptors of nervous terminations) which transmit this stimulation to the spinal cord acting consequently on the specific sensitivity ways for pain.

This stimulation causes an analgesic effect, by stimulating the release of the endorphins so that the pain is relieved naturally.

Destination of use:

Clinical purpose: Therapeutic, for pain relief (TENS)

Field of use: Domestic/Outpatient

The device is indicated for pain relief and for muscular trauma recovery, contractures and distortions. You can treat effectively diseases both of painful and inflammatory component, that benefit from the use of the pulsed current.

Package content:

- **Body Trainer** device
- Protective cloth case
- N. 1 battery 1,5 V AA
- N. 1 Instruction booklet

Warnings:

- Check the position and meaning of all the labels on the device;
- Do not use the device if you see damages and in particular to the casing and to the treatment head: please contact the manufacturer as reported in "Assistance" paragraph;
- Avoid the use of the device to people not properly informed by reading the instruction booklet;
- Avoid the use in wet ambients;
- During the treatment, we recommend to the user not to wear metal objects;
- It is forbidden to place the electrodes in such way that the flow of current crosses the cardiac area;
- It is forbidden the use of the device with the electrodes placed on or near lesions or lacerations of the skin;
- It is forbidden to place the electrodes on the carotid sinus (carotid), genitals;
- It is forbidden to place the electrodes near the eyes and do not invest the eyeball with the output current (an electrode diametrically opposite to the other respect to the eye); keep a minimum distance of 3 cm. from the eyeball;
- Do not use the device on patients with an electronic device implanted without consulting the doctor;
- The simultaneous connection of a patient to a HF surgical device can lead to burns in correspondence to the electrodes of the stimulator and the stimulator can be damaged;
- Operating close (p.e. 1 mt) to a shortwave or microwave device for therapy can produce instability in the output of the stimulator;

The manufacturer is responsible for the performance, reliability and safety of the device only if:

- eventual additions, modifications and/or repairs are performed by authorized personnel;
- the device is used in strict accordance with the instruction of use contained in this manual;

Electromagnetic interference

The device does not generate and does not receive interferences from other equipment. Anyway, it is appropriate to use the device keeping the applicator at least 3 meters away from televisions, monitors, mobile phones or any other electronic devices.

CONTRA-INDICATIONS

Pregnant patients, tuberculosis, juvenile diabetes, viral diseases (acute), mycosis, dermatitis, patient with heart diseases, serious arrhythmias or patient with pacemakers, children, magnetizable denture wearers, acute infections, open wounds, epilepsy (unless different prescriptions).

There aren't any known significant side effects. In some cases on people who are particularly sensitive redness occurred after the treatment corresponding to the application head: the redness usually disappears within few minutes after the treatment.

If the redness persists consult a doctor. In rare case the evening stimulation causes on someone a delay in falling asleep. In such cases avoid the evening treatment.

Body Trainer operates on a battery (1 x AA 1,5 V) it has adjustable intensity from 1 to 10, (to adapt it to your own body) wireless, it does not require the use of electrodes, gel or control unit.

It can be easily used by right or left handed and the ergonomic design has been studied to be easily grabbed at any age and you can reach any area of the body.

You can use it directly on your skin or over light clothes wherever you are, home, car (not driving), boat, on vacation.

It is a valuable help for people that practise sport, before, during and after any sport activity as it relaxes muscles and relieves muscle soreness.

HOW TO USE (for external use only)

We recommend that you always adjust the intensity for every application:

after having properly placed the head on the area to be treated, start from position 1 and gradually adjust the intensity tolerable on your body. In this way you avoid surprise reaction using **Body Trainer** with a too high intensity.

After you have adjusted the intensity, properly place the head on the area to be treated and start doing a "back and forth" movement moving slowly (as you would use an electric shaver). If you use **Body Trainer** over light clothes use it with higher intensity.

Repeat the treatment several times a day with application of at least 10/15 minutes at a time.

WHAT KIND OF PAIN CAN BODY TRAINER RELIEVE?

We can say that **Body Trainer** can be used for various types of pain as:

- Muscle pain
- Back pain
- Shoulder pain
- Upper and lower limbs pain
- Headache and migraine
- Tennis elbow
- Sport traumas and sciatica
- Muscles relaxation and against muscular contractions

This method is used since many years, successfully, in hospitals and pain clinics worldwide.

WHEN NOT TO USE BODY TRAINER:

- If you have a pacemaker
- In the first 16 weeks of pregnancy
- On and around the eyes
- On damp or wet skin
- Near a metal prosthesis
- Near oxygen or flammable liquids
- If you suffer from epilepsy
- On open wounds
- On aortic arteries, on the muscles of the larynx and pharynx
- Over the heart
- While driving any vehicle

CAUTIONS

Read this instruction manual before use.

Keep out of reach of children.

These guidelines do not intend give any medical advice, when in doubt, before using the product and if you are using medications such as steroids, warfarin, etc, you need to consult your doctor.

The use of **Body Trainer** is forbidden to children under the age of 12 years without adult supervision.

Keep out of reach of children; small parts could be removed and swallowed.

INSTRUCTIONS FOR CLEANING

Use a damp cloth. For the treatment head use a damp cloth or non-aggressive detergents.

IMPORTANT




If you suffer from pains or other symptoms that you do not know the origin of, consult your doctor.

Body Trainer can relieve a pain, but in no case can avoid a doctor's diagnosis or replace a treatment without consulting a doctor.

BATTERY CHANGE (1 x AA 1,5 V)

Do not forget to change the battery in case the device does not work.

- 1 - Slide the battery cover as in illustration A to open it.
- 2 - Remove the old battery and put the new battery AA (1,5 V) as in illustration B.
- 3 - Slide the battery cover as in illustration C to close.

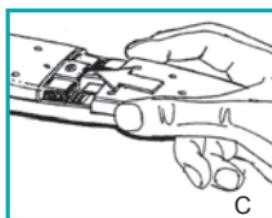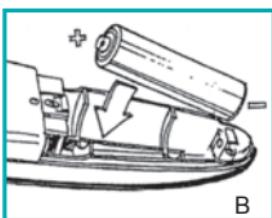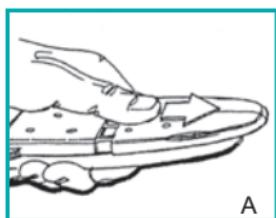Body Trainer has **2 years warranty** (with replacement) against any manufacturing defect, with regular and proper use of the device.

The warranty decays in case of tampering and in case of intervention on the device by personnel non authorized by the manufacturer or by the authorized retailer.

CAUTIONS FOR TRANSPORT

There are no specific care to use during the transport because **Body Trainer** is a portable device.

We recommend to put **Body Trainer** after every use in the supplied case.

CAUTIONS FOR STORAGE

The equipment is protected up to the following environmental conditions:

Without the supplied package	
ambient temperature relative humidity pressure	from +5 to + 40 °C from 30 to 80% from 500 to 1060 hPa
With the supplied package	
ambient temperature relative humidity pressure	from -10 to + 50 °C from 10 to 90% from 500 to 1060 hPa

INFORMATION FOR DISPOSAL

The product is subject to the WEEE regulation (presence of the symbol on the label) on recycling: for the disposal of the product, please use the specific equipped areas for the collection of electronic material or contact the manufacturer.

MAINTENANCE AND PROBLEM SOLVING

If the equipment is used according to what indicated in this manual, any particular routine maintenance is required.

If you encounter malfunctions or problems in the use of **Body Trainer**, check the following:

- **It does not turn on.** Check the charging status of the battery and eventually replace it (see chapter "Battery change"). If the problem continues contact the manufacturer.
- **It does not transmit electrical impulses.** Check that you have properly turned the regulation wheel and check that the led light is on. If the problem continues contact the manufacturer.
- **Body Trainer turns off during operation.** We recommend to replace the battery and start again the treatment. If the problem continues contact the manufacturer.

ASSISTANCE

The manufacturer is the only exclusivist for technical intervention on the equipment.

For any technical assistance please contact:

AXEL S.r.l.

Via F. Matteucci, 11 - 50041 Calenzano (FI) Italy


Tel. 055 8827245 - Fax 055 8826338

Eventual technical documentation concerning repairable parts can be provided, but only upon company authorization and only after giving adequate instructions to the intervention personnel.

SYMBOLS

Manufacturer: AXEL s.r.l. - Via F. Matteucci, 11 - 50041 Calenzano (FI)

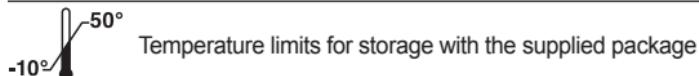Temperature limits for storage with the supplied package

Follow the instruction for use

Warning

WEEE directive warning:

old electrical appliances are valuable materials and should not be disposed of with domestic refuse! Customers are urged to help protect the environment and natural resources by taking the appliance to a special waste collection centre when it is no longer used.

Compliance Medical Devices Directive - CEE 93/42 - CE 2007/47

Type BF applied part

Use in a dry place
and away from water

TECHNICAL SPECIFICATIONS

Description: portable cutaneous electrostimulator

Product dimensions: 220 x 50 mm

Classification applied part: Type BF

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© Copyright, 2011, AXEL s.r.l. Tutti i diritti riservati. Copie o altre riproduzioni di questo documento sono proibite, previa autorizzazione scritta della AXEL s.r.l. Italy.
© Copyright, 2011, AXEL s.r.l. All rights reserved. Copies or reproductions of this document are forbidden unless specifically authorised in writing by AXEL s.r.l. Italy.